


Agilent Bio SEC-5 Guard Column - New in Box

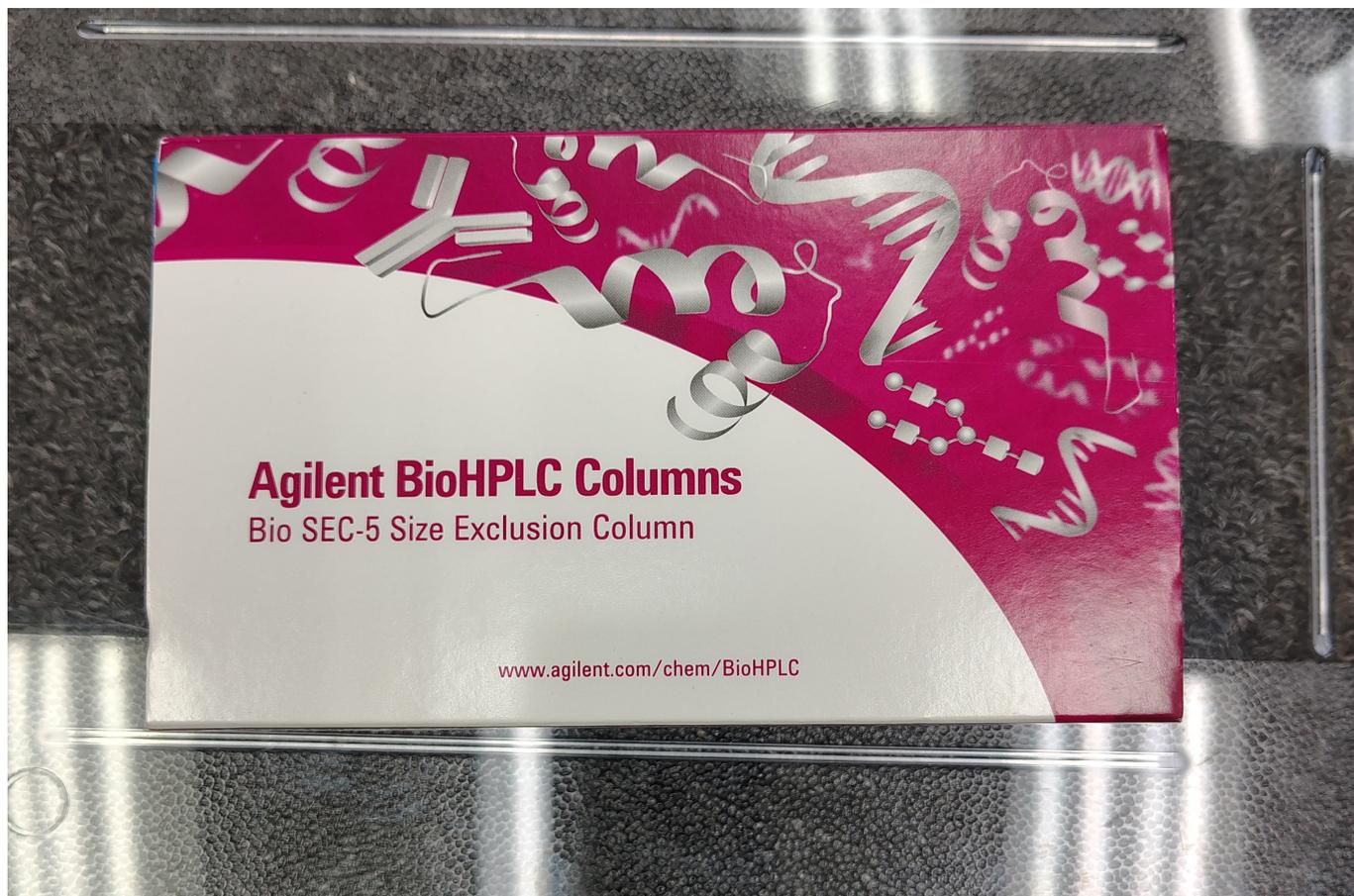

Particle size: 5 μm ; Pore size: 300 \AA

Rating: Not Rated Yet

Price

Sales price \$250.00

[Ask a question about this product](#)

Manufacturer [Agilent](#)

Description

The Agilent Bio SEC-5 Guard column is a size-exclusion chromatography (SEC) guard column designed to protect and extend the lifetime of a main analytical SEC column in biomolecule applications. It uses 5 μm spherical, fully-porous silica particles coated with a neutral hydrophilic polymeric layer, enabling stable and reproducible separations of proteins, peptides, and other macromolecules.

Specifications

Particle size: 5 µm

Pore size: 300 Å

Column dimensions: 4.6 mm × 50 mm length

pH operating range: approx. pH 2.0 to 8.5.

Comes with our 90-day warranty.

Biotech Equipment Sales, Inc.: based in South San Francisco, the Biotech capital of the world, since 2009. We are dedicated to bringing high-quality pre-owned laboratory instrumentation to our scientific community.

Biotech Equipment Sales, Inc.
226 Miller Ave.
South San Francisco, CA 94080
650-871-5707 office
650-276-7487 fax
sales@biotechequipmentsales.com