


## Agilent BioTek 405 TS Touch Microplate Washer 405TSUV

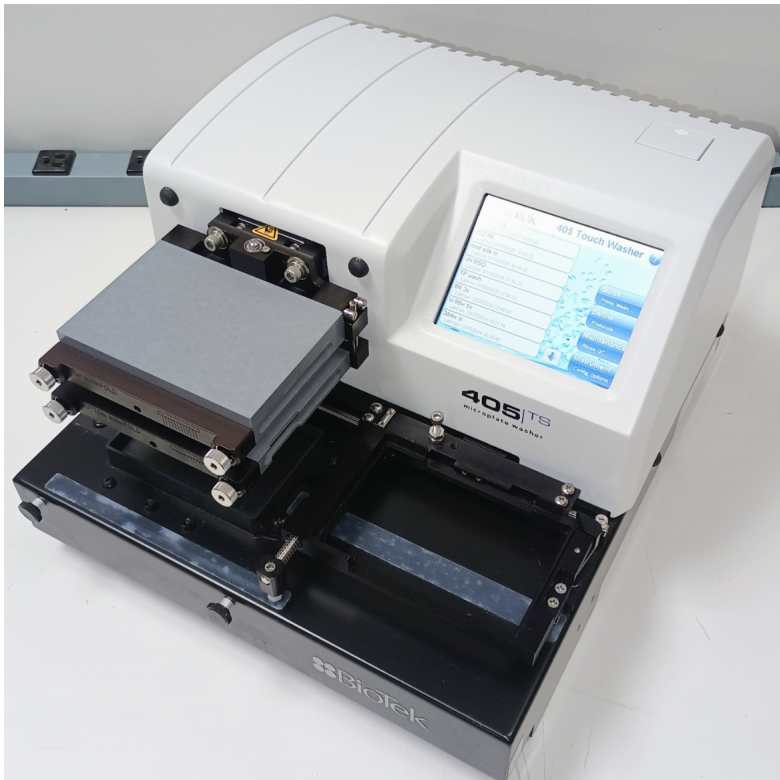

For both 96- and 384-well microplate washing.

Rating: Not Rated Yet

**Price**

Sales price \$9,500.00

[Ask a question about this product](#)

Manufacturer [Agilent](#)

**Description**

The Agilent BioTek 405TS uses a high resolution touch screen interface for all programmable operations. Features programmable dispense volumes, flow, and aspiration rates and provides for a wide range of washing capabilities. Supports wash, prime, dispense, and aspirate programs. Can be used to wash all microplate-based assays, including ELISA, fluorescence, chemiluminescence, RNA, DNA probes, and cellular assays.

**Features:**

"U" Version is for both 96- and 384-well microplate washing

"V" Version allows for up to 4 buffer bottles to be used with its built-in valve assembly

100-240V, 50/60Hz, 9A. What is

Included:

405TS 405TSUV microplate washer

Standard Thomas/Gast/Welch piston vacuum pump (Upgrade to Vacuubrand pump available)

Two 4-Liter waste bottles with vacuum tubing (Upgrade to 20-Liter primary waste bottle available)

One 4-Liter supply bottle with dispense tubing (Upgrade to full set of 4 supply bottles available)\*Biostack plate stacking systems are also available

*Comes with our 90-day warranty.*

Biotech Equipment Sales, Inc.: based in South San Francisco, the Biotech capital of the world, since 2009. We are dedicated to bringing high-quality pre-owned laboratory instrumentation to our scientific community.

**Biotech Equipment Sales, Inc.**  
**226 Miller Ave.**  
**South San Francisco, CA 94080**  
**650-871-5707 office**  
**650-276-7487 fax**  
[sales@biotechequipmentsales.com](mailto:sales@biotechequipmentsales.com)