


Cole Parmer MasterFlex Peristaltic Pump

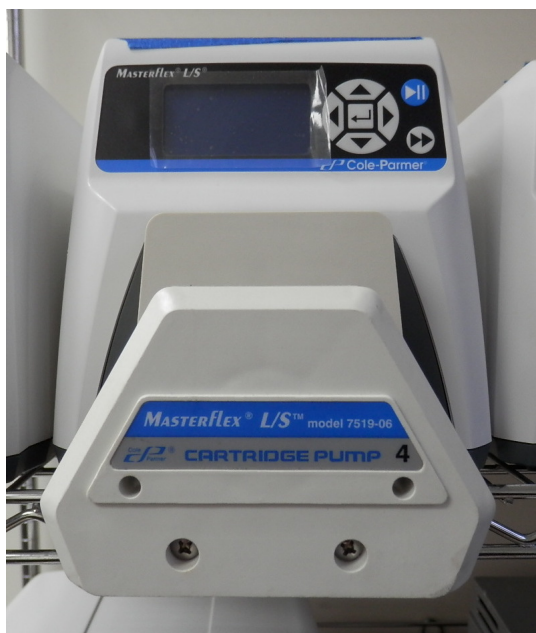

Four operating modes: continuous run, timed dispense, volume dispense, and copy dispense, and with remote capabilities.

Rating: Not Rated Yet

Price

Sales price \$1,500.00

[Ask a question about this product](#)

Manufacturer [Cole-Parmer](#)

Description

The Cole Parmer Masterflex L/S is a Digital Drive Peristaltic Pump with digital control and a wide speed range. This Cole Parmer Pump has a graphical interface that allows users to select and program four operating modes: continuous run, timed dispense, volume dispense, and copy dispense, and with remote capabilities, this digital drive is an ideal choice for automated process applications. The dispense interval also allows users to set a delay, and set on-off cycles, and features an anti-drip functionality to ensure dispensing accuracy.

The Pump runs on standard 120V Power and can be used with several different pump heads.

Specifications:

Control Type: Digital Variable Speed

Dispense by Volume: 0.001 to 9,999 mL, L, or gal; 1 to 9,999 dispense cycles

Time interval of on/off cycles (1 sec to 99 hr, 59 min, 59 sec per cycle)

Min RPM (rpm): 0.02 (7523-90); 0.1 (7523-80)

Max RPM (rpm): 100 (7523-90); 600 (7523-80)

Speed Control $\pm 0.1\%$ (0.1 rpm at 600 rpm): Variable

Reversible Motor: Yes

Digital Dispensing: Yes

Peristaltic Pumps: Cole Parmer MasterFlex Peristaltic Pump

No. of Pump Heads Accepted: 2

Type of Pump Heads Accepted: Standard, Easy-Load, Easy-Load II, Easy-Load 3, High Performance, Cartridge, Multichannel, PTFE Tubing, PTFE Diaphragm

Drive IP rating: IP33

Motor: hp1/10

Remote I/O: 0 to 20 mA, 4 to 20 mA, 0 to 10 V, start/stop, direction, and prime via contact closure

Dimensions: 9x15x9" (W x D x H)

Power: 120V 50/60Hz 2.2A

Weight: 15lbs

Comes with our 90-day warranty.

Biotech Equipment Sales, Inc.: based in South San Francisco, the Biotech capital of the world, since 2009. We are dedicated to bringing high-quality pre-owned laboratory instrumentation to our scientific community.

Biotech Equipment Sales, Inc.
226 Miller Ave.
South San Francisco, CA 94080
650-871-5707 office
650-276-7487 fax
sales@biotechequipmentsales.com