


Sartorius MSA 225S100DI

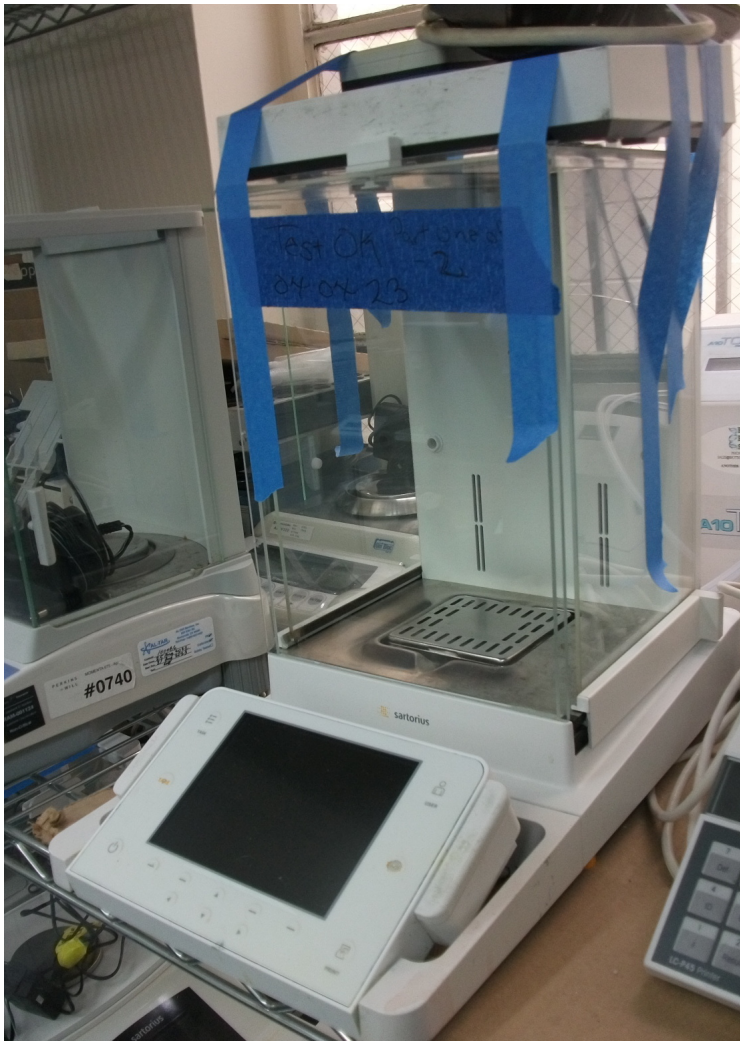

Capacity Resolution: 220 g x 0.01 mg; Reproducibility: (\pm) 0...60 g: 0.015 mg | 60...220 g: 0.025 mg

Rating: Not Rated Yet

Price

Sales price \$4,750.00

[Ask a question about this product](#)

Manufacturer [Sartorius](#)

Description

Balances: Sartorius MSA 225S100DI

Sartorius Cubis Semi-Micro Balance with MSA Display and DI Draft Shield

Capacity Resolution: 220 g x 0.01 mg

Reproducibility: (\pm) 0...60 g: 0.015 mg | 60...220 g: 0.025 mg

Linearity: (\pm) 0.1 mg

External Calibration Point: 200g

Data Interface: Ethernet

Chamber Dimensions: 154mm(Depth) x 205mm(Width) x 250mm(Height)

Features:

Auto door feature, open and close motorized doors with keypad, touchless switch, foot switch or commands from PC.

Includes integrated ionizer to eliminate static from inside the weighing chamber, recommended when dealing with plastics or powders. Unlike other ionizers this does not use a fan or radioactive material.

This balance includes motorized leveling, the balance senses an auto level condition and self corrects with motorized feet

Comes with our 90-day warranty.

Biotech Equipment Sales, Inc., has been based in South San Francisco, the Biotech capital of the world, since 2009. We are dedicated to bringing high-quality pre-owned laboratory instrumentation to our scientific community.

Biotech Equipment Sales, Inc.
226 Miller Ave.
South San Francisco, CA 94080
650-871-5707 office
650-276-7487 fax
sales@biotechequipmentsales.com