


Thermo Forma Steri-Cycle CO2 Incubator

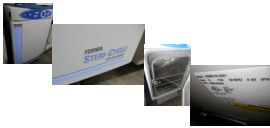

6.5 cu ft air jacketed single chamber CO2 incubator.

Rating: Not Rated Yet

Price

Sales price \$3,250.00

[Ask a question about this product](#)

Manufacturer [Thermo Forma](#)

Description

This 6.5 cu ft air jacketed single chamber CO2 incubator uses HEPA filtration to cycle the air completely every 60 seconds for cleanroom Class 100 conditions. The incubator also features a high heat sterilization cycle, running at 140°C over a 14-hr. cycle for easy cleaning and sterilizing, with a check prompt to prevent accidental activation. It also features fan-assisted airflow with an IR CO2 sensor for a uniform growth environment. The interior is a mirror finish stainless steel for easy cleaning, with a stainless humidity pan and removable stainless shelving that's easily placed without the need for extra tools.

The Cell Culture incubator has a tempered glass door to help maintain the environment while allowing visual monitoring. The microprocessor control gives you a centralized command center up front, with a bright LED display, and allows users to choose their desired temperature range

Incubators: Thermo Forma Steri-Cycle CO2 Incubator

between 5°C above ambient to 50°C, humidity, and CO2 concentration, as well as set alarm parameters.

Comes with our 90-day warranty.

Biotech Equipment Sales, Inc.: based in South San Francisco, the Biotech capital of the world, since 2009. We are dedicated to bringing high-quality pre-owned laboratory instrumentation to our scientific community.

Biotech Equipment Sales, Inc.
226 Miller Ave.
South San Francisco, CA 94080
650-871-5707 office
650-276-7487 fax
sales@biotechequipmentsales.com