


## Perkin Elmer Victor X4 Multilabel Microplate Reader for TRF

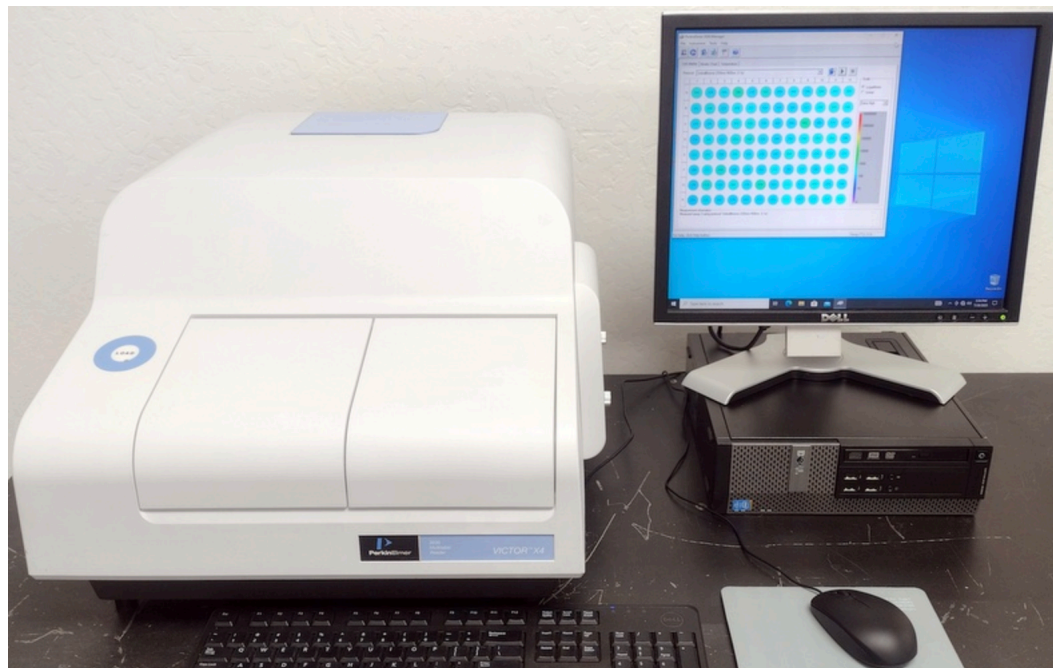

Comes with a Windows 10 computer.

Rating: Not Rated Yet

**Price**

Sales price \$13,990.00

[Ask a question about this product](#)

Manufacturer [Perkin Elmer](#)

**Description**

Version 2030-0040 Time-Resolved Fluorometry model, capable of fluorometry (top and bottom reading), luminometry (flash and glow), photometry (UV and Visible), time-resolved fluorometry, and LANCE, installed hardware options include temperature control, adjustable beam size, counting from below, and UV absorbance, comes with a Windows 10 computer with 2030 Workstation Manager v.4 installed, monitor, keyboard and mouse, 110-120V/220-240V, 50/60Hz.

Standard excitation and emission default filter sets (A) are included which covers absorbance at 405nm, 450nm, and 490nm and fluorescence at 355/460nm (umbelliferone) and 485/535nm (fluorescein). Custom filter/s can be modified depending on assay requirements.

*Comes with our 90-day warranty.*

Biotech Equipment Sales, Inc.: based in South San Francisco, the Biotech capital of the world, since 2009. We are dedicated to bringing high-quality pre-owned laboratory instrumentation to our scientific community.

**Biotech Equipment Sales, Inc.**  
226 Miller Ave.  
South San Francisco, CA 94080  
650-871-5707 office  
650-276-7487 fax  
[sales@biotechequipmentsales.com](mailto:sales@biotechequipmentsales.com)