


VWR Cool Thermal Shake Touch Thermoshaker

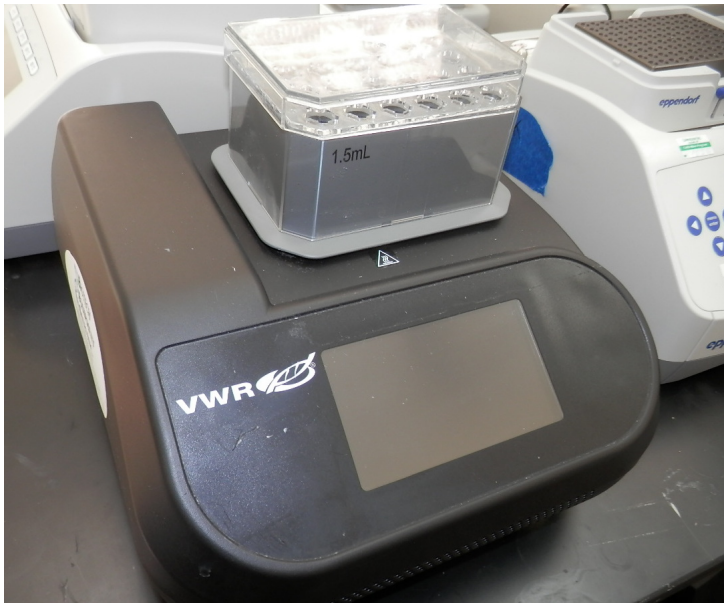

Heats, shakes, cools. With intuitive 4.3" LCD Touch Screen

Rating: Not Rated Yet

Price

Sales price \$1,425.00

[Ask a question about this product](#)

Manufacturer [VWR](#)

Description

Designed for applications that require consistent and precise results.

Intuitive 4.3" LCD Touch Screen simultaneously displays all settings for complete control and monitoring

Rapid temperature cycling and high-speed shaking ability result in increased throughput

Internal memory stores five separate 5-step programs, and operational data can be exported via USB

Interchangeable blocks are ideal for processing micro tubes, conical tubes, cryotubes, and microplates

With heating, cooling, and shaking capabilities, the Cooling Thermal-Shake Touch uses interchangeable blocks to accommodate a wide variety of tubes and microplates. The low-profile design minimizes the unit's footprint on the bench. The LCD touch screen enables faster input of temperature, speed, and time, which can all be viewed at once. Display features on-screen help topics with operational tips and touch screen is compatible with rubber gloves used in labs. USB port can transfer information to a flash drive for data logging, program storage, and software updates. Program control capabilities allow user programmable operation for automated use plus storage of five separate 5-step programs. Adjustable temperature ramp rate feature separately defines the heating and cooling rate in increments of 0.5 degrees/min. For maximum temperature accuracy, the single point calibration procedure allows the user to calibrate at 6 distinctly defined temperatures. Though the unit can process an array of sample sizes in various VWR Blocks (sold separately), the Cooling Thermal-Shake Touch conveniently and effectively operates with non-VWR Blocks. The unit is equipped with a pulse mode feature for quick vortex applications.

Comes with our 90-day warranty.

Biotech Equipment Sales, Inc.: based in South San Francisco, the Biotech capital of the world, since 2009. We are dedicated to bringing high-quality pre-owned laboratory instrumentation to our scientific community.

Biotech Equipment Sales, Inc.
226 Miller Ave.
South San Francisco, CA 94080
650-871-5707 office
650-276-7487 fax
sales@biotechequipmentsales.com